

Freely movable DUCTLESS FUME HOOD

• Carbon filter (option)
HEPA filter (option)
Pre-filter

Anemometer

It is available to check easily air velocity inside the chamber.

Filter test hole

Saturation state of the filter can be checked by use of the gas detection Kit.

LCD controller

Product condition can be easily observed, and operation of blower/lamp etc. is simple.

Three sides visibility with clear acrylic / PMMA

Securing of the optimum view is available as the transparent acrylic is used for the three sides to get an easy observation of the work surface.

Utility hole


Incoming of electricity, gas line into the work space is free as the utility hole is equipped.

Water tap / Cupsink (Option)

Water tap and cupsink can be additionally selected for easy washing during the experiment.

Leveling caster

It can be easily moved to any desired places as the caster of stopper function is attached.

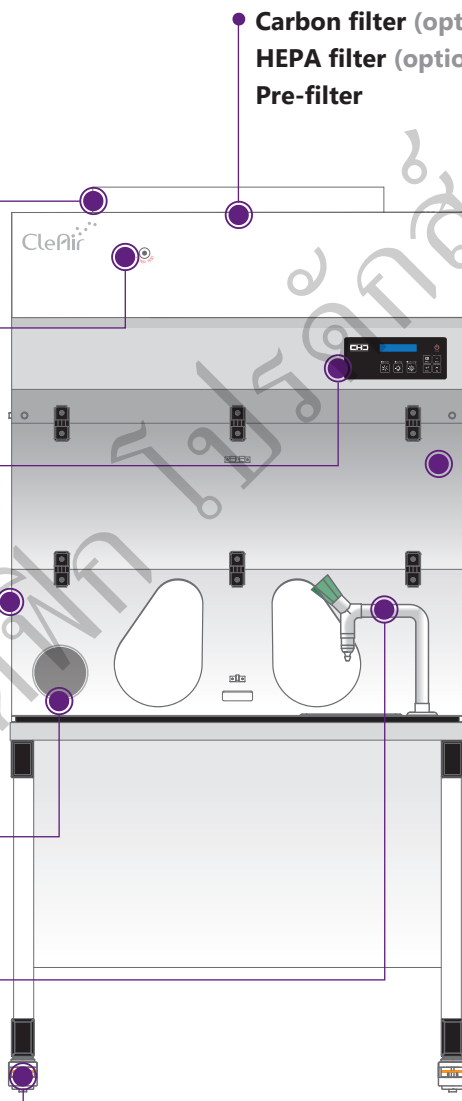

• Carbon filter type

Type	Applications
Organic filter	Fume(gas, vapor) from organic solvents
Acid gas filter	Fume(gas, vapor) from acids



CLE-051, CLE-051F

Ductless Fume Hood



CLE-051F-03



Airflow

- Room air
- Fume (gas, vapor)
- Clean air



General specification		Model No.	CLE-051		CLE-051F	
			03	04	03	04
Dimension	Exterior	W/mm	900	1200	900	1200
		D/mm	660	660	660	660
		H/mm	1100	1100	1970	1970
	Interior	'W/mm	884	1184	884	1184
		'D/mm	625	625	625	625
		'H/mm	730	730	730	730

※ CLE-051 excludes the support stand.